
ORBI-Sept Wet Wipes

Nachfüllpackung | Refill

Soft M - Duftnote neutral

- de Desinfektionstücher Nachfüllpackung
- en Disinfection Wipes Refill
- da Desinfektionsservietter – genopfyldningspakke
- sv Desinfektionsservetter, påfyllningsförpackning
- no Desinfeksjonsservietter, påfyllingspakke
- fi Desinfiointipyyhkeiden täyttöpakkaus
- pl Chusteczki dezynfekujące – opakowanie uzupełniające
- nl Navulverpakking met desinfectiedoekjes
- cs Náhradní balení dezinfekčních ubrousků
- fr Recharge de lingettes désinfectantes

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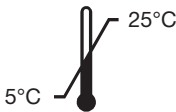Types of application:

Ready-to-use wipes for rapid prophylactic disinfection of alcohol-resistant surfaces of sensitive medical devices. Intended for use in the medical sector ((dental) medical practices, hospitals, nursing homes, care and rescue services) as well as in the food sector, industry and public facilities.

Disinfection of non-invasive medical devices:

Preparation for prophylactic wipe disinfection of alcohol-resistant medical devices such as dental treatment units, hospital beds, blood pressure monitors, venous tubes, handpieces and contra-angles etc. To disinfect, wipe the surface of the medical device with the soaked cloth.

Instructions for use:

For labelling, remove the stick & lift label from the bag and stick it on the dispensing container. Open the bag at the top and place it in the dispensing container. Pull on the first cloth of the fleece roll and thread it through the star-shaped lid. Remove the following cloths. After removing the cloth, press the lid firmly onto the container and close it tightly. Once the cloths have been used up completely, the dispensing container must be cleaned and disinfected in accordance with the reprocessing instructions. Reprocessing should be carried out with a suitable, preferably alcohol-based disinfectant wipe. All surfaces to be disinfected must be completely wetted with the disinfectant solution. Allow the product to take effect according to the specified contact times. Only use on alcohol-compatible materials (e.g. not suitable for acrylic glass). Use is intended for professional users only. The wipes can be used for a maximum of 28 days after the bag has been opened for the first time. The material compatibility of new or unknown materials should be tested on a non-visible area.


Contact times for the prophylactic disinfection of medical devices and surfaces in the medical sector:

Bacteria and Candida albicans according to EN 13727, EN 13624, EN 16615 and VAH	Conc. 1 min.
under dirty conditions:	Conc. 1 min.
Tb (tuberculosis) inactivation (M. terrae, EN 14348, dirty conditions):	Conc. 1 min.
„Limited virucidal“ (enveloped viruses such as e.g. HBV, HIV, HCV) according to EN 14476	
and DVV/2012, dirty conditions:	Conc. 1 min.
„Limited virucidal PLUS“ (enveloped viruses such as HBV, HIV, HCV) and non-enveloped Noro,	
Adeno and Rota viruses according to EN 14476 (dirty conditions):	Conc. 1 min.

Safety instructions:

H 228 Flammable solid. P210 Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. P280 Wear protective gloves / eye protection. P404 Store in a closed container.

All serious incidents related to the device shall be reported by the user/patient to the manufacturer and the competent authority of the Member State where the user/patient is established.

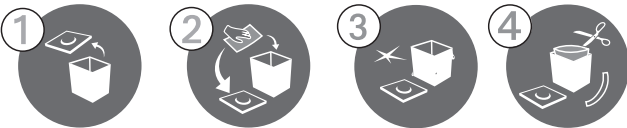Preparing the dispenser box:

1. Remove the lid from the container. 2. Wipe the lid, inside of the container, wipe dispenser and reservoir with a soaked disinfectant wipe (observe the exposure time of the disinfectant). 3. Allow the container and lid to dry. 4. Cut open the filled bag at the top and insert it into the container. The disinfectant wipes are ready for use.

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